



TrueBond

TrueBond

Multipurpose Veterinary Centrifuge

Part # 00-376-009-005

Using the TrueBond multipurpose veterinary centrifuge, vet laboratories can now spin microhematocrit tubes, fecal, blood, urine, and more in one centrifuge.

Spin Any Tube

TrueBond is designed to spin all your tubes, from microhematocrit to urine and everything in between. Includes all necessary adapters.

Easy to Use

Minimize training time with the user-friendly interface. The fully digital display makes it easy to customize cycle parameters on the fly, or set and select one of the 10 programmable cycles.

Easy to Monitor

Process samples more efficiently with TrueBond's LED lid lighting. Once the unit is started, the display shows remaining spin time.

Safe and Durable

Enjoy years of safe, worry-free operation with the brushless motor, imbalance protection, and a clear, shatterproof lid, which locks until the centrifuge has stopped completely.

A Great Value


Each TrueBond centrifuge is ready to run with fixed angle rotor, two sets of tube holders, microhematocrit carriers, microtainer inserts, and 2-year warranty included in the affordable price.



Applications	Max Speed	Max Capacity	Included	Warranty
Microhematocrit, Urine, Blood, Fecal	1,850 xg Max RCF 3,900 Max RPM	6x fecal, urine, and blood tubes (.5-15 mL) 8x microhematocrits	Rotor (Fixed Angle) 6x .5-10 mL Tube Holders 6x 15 mL Tube Holders 2x Microhematocrit Carriers	2 years Rotor: lifetime warranty

Designed, built, and supported in the USA

Free 30-Day Evaluation (U.S. Only)

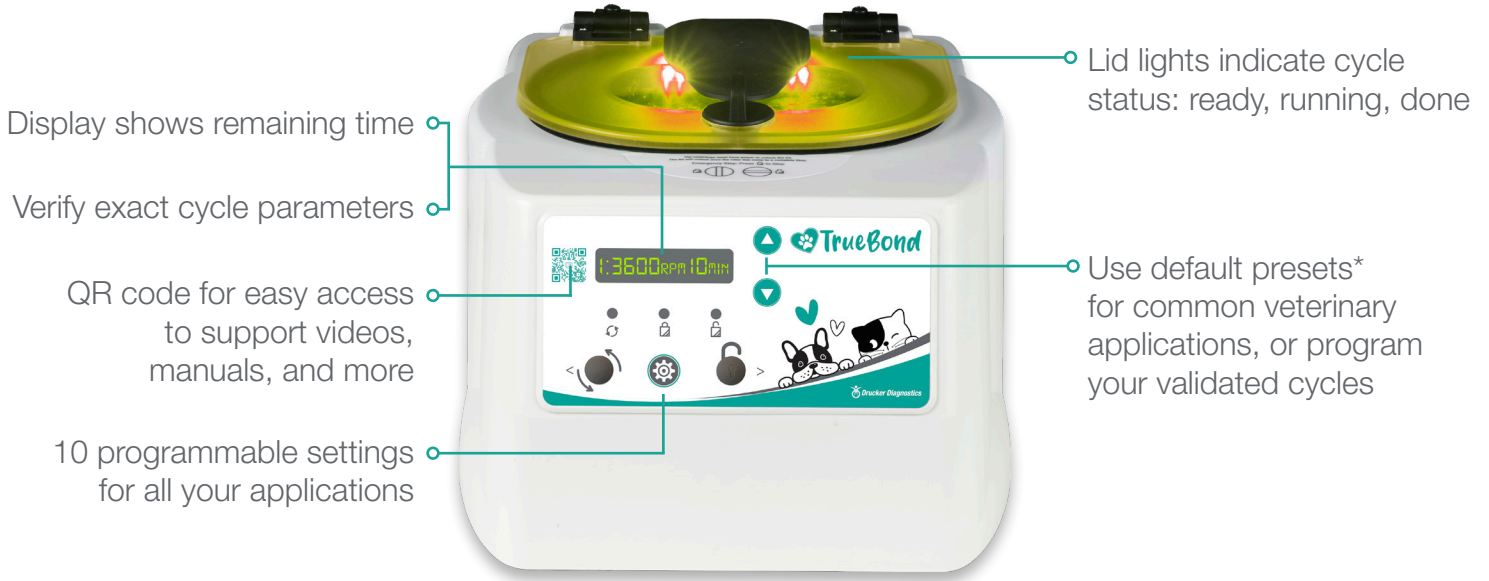

Standard Preset Cycles		Blue Tube Holder (75-100 mm/3-10 mL tubes)	Orange Tube Holder (125 mm/15 mL tubes)
1 (Chemistry)	10 min @ 3,900 RPM	1,600 xg	N/A
2 (Urine Sediment)	5 min @ 2,000 RPM	N/A	500 xg
3 (Veterinary Fecal)	5 min @ 1,600 RPM	300 xg	N/A
4 (Microhematocrit)	6 min @ 3,800 RPM	1,550 xg	N/A

Included With TrueBond Purchase	
Centrifuge	00-376-009-005
75-100 mm Blue Tube Holders	7713079K (pack of 6)
125 mm Orange Tube Holders	7713044K (pack of 6)
Microhematocrit Carriers	02-004-0-0027K (pack of 2)
1-1.5 mL Orange Inserts	7713065K (pack of 6)
0.5-1 mL Purple Inserts	7713068K (pack of 6)

Specifications	
Sound Level	64 dBA
Max RPM/G-Force	3900 RPM
Max G-Force (Orange Tubes)	1600 xg
Max G-Force (Blue Tubes)	1850 xg
Dimensions (W x D x H)	12 x 14 x 9 in 30 x 36 x 23 cm
Weight	12 lbs (5.4 kg)
Power	95 to 253 V – 50/60 Hz



druckerdiagnostics.com
sales@druckerdiagnostics.com

+1-866-265-1486 (U.S. only)
+1-814-692-7661

Protected by U.S. Patents #6.811.531. #D718,463S #D734,489S Other Patents Pending.